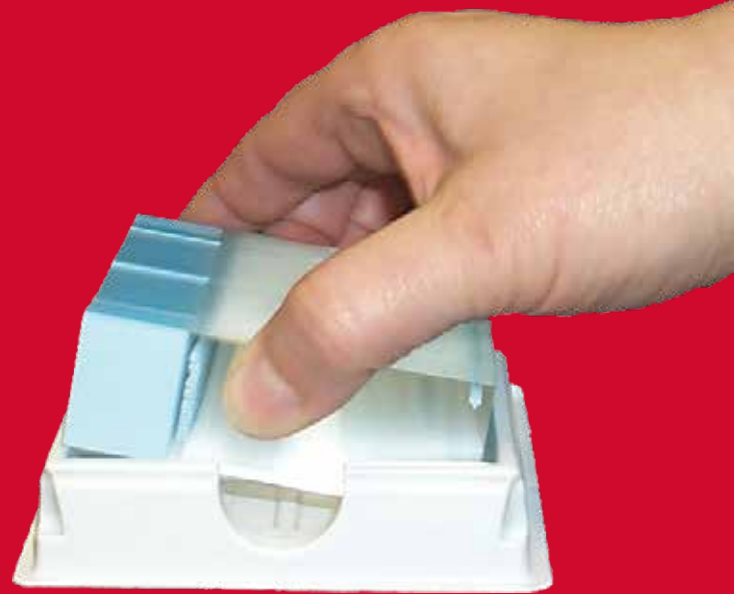


Histology Innovation for a NEW Generation

Staining & Coverslipping

SHURMark® Microscope Slides, Cover Glass & Accessories



**One Source For All Your Laboratory Needs.
Protecting Every Life Story In Your Lab.**

Each tissue specimen your lab processes has a life story behind it. Your mission is to ensure your lab has the best instruments and resources in order to deliver optimal clinical outcomes and the best possible care. Our comprehensive microscope slides & cover glass product line includes everything your lab needs for increased operational efficiency, advanced specimen safety and highly accurate results.



High-Quality Clarity Delivers Optimal Light Transmission

SHUR**Mark**®s microscope slides are manufactured from high quality clear borosilicate float glass (refractive index of $n_e=1.51$), hydrolytic Class I, to exacting specifications in a modern high-tech manufacturing facility. Automated ultrasonic cleaning tanks, ultra-pure water and infrared drying ovens are utilized for maximum cleanliness & dryness. Final product is then hermetically sealed & desiccant added to maintain freshness.

Select the ideal slide type for your application based on your method of applying specimen identification information; i.e.- pencils, pens, automated ink-jet or scribing labelers.

Variety of Slide Types to Fit Every Labs Application

- Colors:** 6 popular chemical resist colors
- Coating:** Available in either uncoated or coated for firm specimen attachment
- Slide Corner:** Select either a 90° 45° or Beveled corner
- Application Method:** Apply specimen identification information via pen, pencil, automated ink-jet or scribing labelers.

Minimized Jamming in Automated Labelers Provides Operator Safety

All sizes have excellent flatness, chemical resistivity and optimal light transmission values. Each edge is precision ground to a smooth finish for operator safety.

Moisture Resistant Packaging Minimizes Sticking

All microscope slides are (1 x 3) (25 x 75)(in)(mm); packaged 72/pk, 20pk/cs in a moisture resistant packaging to minimize sticking.

No. 1

Cat. #	Description 1oz/ea; 10ea/pk	Thickness (mm)	# of pcs/ea
SM-1818-1	(18x18)mm	0.13 - 0.17	233
SM-2020-1	(20x20)mm	0.13 - 0.17	190
SM-2222-1	(22x22)mm	0.13 - 0.17	156
SM-2230-1	(22x30)mm	0.13 - 0.17	114
SM-2232-1	(22x32)mm	0.13 - 0.17	103
SM-2240-1	(22x40)mm	0.13 - 0.17	87
SM-2250-1	(22x50)mm	0.13 - 0.17	70
SM-2260-1	(22x60)mm	0.13 - 0.17	58
SM-2424-1	(24x24)mm	0.13 - 0.17	130
SM-2430-1	(24x30)mm	0.13 - 0.17	105
SM-2432-1	(24x32)mm	0.13 - 0.17	100
SM-2440-1	(24x40)mm	0.13 - 0.17	80
SM-2450-1	(24x50)mm	0.13 - 0.17	65
SM-2455-1	(24x55)mm	0.13 - 0.17	58
SM-2460-1	(24x60)mm	0.13 - 0.17	54
SM-2530-1	(25x30)mm	0.13 - 0.17	104
SM-2540-1	(25x40)mm	0.13 - 0.17	79
SM-2550-1	(25x50)mm	0.13 - 0.17	63

No. 1.5

Cat. #	Description 1oz/ea; 10ea/ pk	Thickness (mm)	# of pcs/ea
SM-2222-1A	(22x22)mm	0.16 - 0.19	134
SM-2240-1A	(22x40)mm	0.16 - 0.19	74
SM-2250-1A	(22x50)mm	0.16 - 0.19	60
SM-2430-1A	(24x30)mm	0.16 - 0.19	90
SM-2440-1A	(24x40)mm	0.16 - 0.19	68
SM-2450-1A	(24x50)mm	0.16 - 0.19	56
SM-2455-1A	(24x55)mm	0.16 - 0.19	50
SM-2530-1A	(25x30)mm	0.16 - 0.19	89

Hermetically
sealed outer
packaging for longer
shelf life in humid
environments

No. 2

Cat. #	Description 1oz/ea; 10ea/pk	Thickness (mm)	# of pcs/ea
SM-1818-2	(18x18)mm	0.19 - 0.23	170
SM-2020-2	(20x20)mm	0.19 - 0.23	138
SM-2222-2	(22x22)mm	0.19 - 0.23	113
SM-2230-2	(22x30)mm	0.19 - 0.23	83
SM-2232-2	(22x32)mm	0.19 - 0.23	75
SM-2240-2	(22x40)mm	0.19 - 0.23	64
SM-2250-2	(22x50)mm	0.19 - 0.23	51
SM-2260-2	(22x60)mm	0.19 - 0.23	42
SM-2424-2	(24x24)mm	0.19 - 0.23	95
SM-2430-2	(24x30)mm	0.19 - 0.23	76
SM-2432-2	(24x32)mm	0.19 - 0.23	73
SM-2440-2	(24x40)mm	0.19 - 0.23	59
SM-2450-2	(24x50)mm	0.19 - 0.23	47
SM-2455-2	(24x55)mm	0.19 - 0.23	42
SM-2460-2	(24x60)mm	0.19 - 0.23	39
SM-2530-2	(25x30)mm	0.19 - 0.23	75
SM-2540-2	(25x40)mm	0.19 - 0.23	58
SM-2550-2	(25x50)mm	0.19 - 0.23	65

Our premium cover glass is made from pure white flat sheets of hydrolytic Class I borosilicate glass (refractive index $n_e=1.525\pm0.0015$) in accordance with Din ISO 8255-1:2003-05. Ultrasonically pre-cleaned, packaged in plastic boxes with moisture absorbing desiccant and then hermetically sealed outer packaging for longer shelf life. No self-fluorescence. Ideal for use in automated coverslipping instruments.












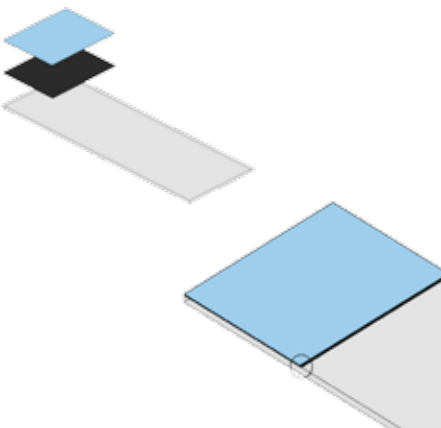
SHURMARK® MICROSCOPE SLIDES FOR SHURMARK® SLIDE ETCHER










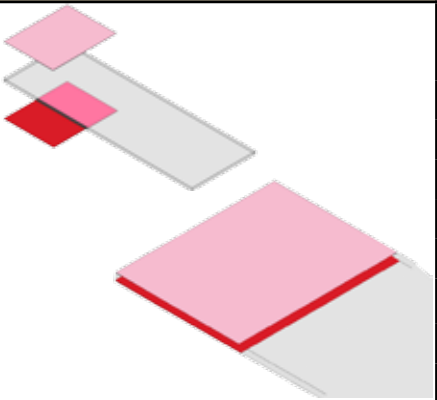
The SHURMark® Slide Etcher utilizes a diamond stylus to permanently etch the paint from the slide. The etched content is not affected by staining, reagents (including Xylene) or physical wear. The SHURMark® Slide Etcher can print on the SHURMark® Contrasting Slide.

Warning: The SHURMark® Slide Etcher will remove paint from a majority of TBS and competitive slides. However, ceramic coated or ink treated slides rapidly increases wear to the diamond stylus and produces inadequate printing.

SHURMark® Contrasting Black Color Under Primary Surface Color for use with Etching Labelers							
Corners		 90° Corners		 45° Clipped Corners		 Beveled Edges	
Colors		UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #
	WHITE	SM-CU-W	SM-CCU-W	SM-CU-45W	SM-CCU-45W	SM-CU-BW	SM-CCU-BW
	BLUE	SM-CU-B	SM-CCU-B	SM-CU-45B	SM-CCU-45B	SM-CU-BB	SM-CCU-BB
	YELLOW	SM-CU-Y	SM-CCU-Y	SM-CU-45Y	SM-CCU-45Y	SM-CU-BY	SM-CCU-BY
	GREEN	SM-CU-G	SM-CCU-G	SM-CU-45G	SM-CCU-45G	SM-CU-BG	SM-CCU-BG
	LILAC	SM-CU-L	SM-CCU-L	SM-CU-45L	SM-CCU-45L	SM-CU-BL	SM-CCU-BL
	PINK	SM-CU-P	SM-CCU-P	SM-CU-45P	SM-CCU-45P	SM-CU-BP	SM-CCU-BP



SHURMark® Contrasting Dark Color on Bottom of Slides for use with Etching Labelers				
Corners			 90° Corners	
Colors			UNCHARGED Cat #	CHARGED Cat #
	Top	Bottom		
	WHITE	BLACK	SM-STT-W	SM-SCTT-W
	BLUE	DK BLUE	SM-STT-B	SM-SCTT-B
	YELLOW	RED	SM-STT-Y	SM-SCTT-Y
	GREEN	DK GREEN	SM-STT-G	SM-SCTT-G
	LILAC	BLACK	SM-STT-L	SM-SCTT-L
	PINK	RED	SM-STT-P	SM-SCTT-P

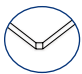

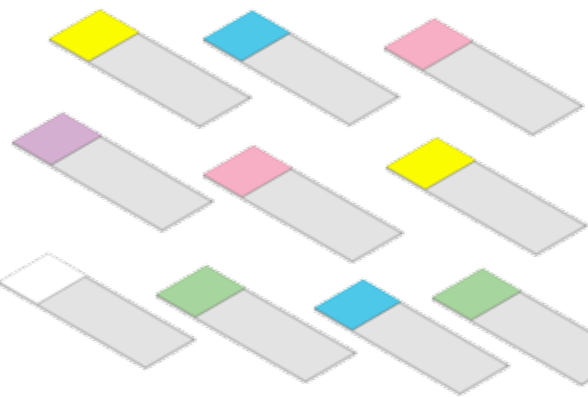










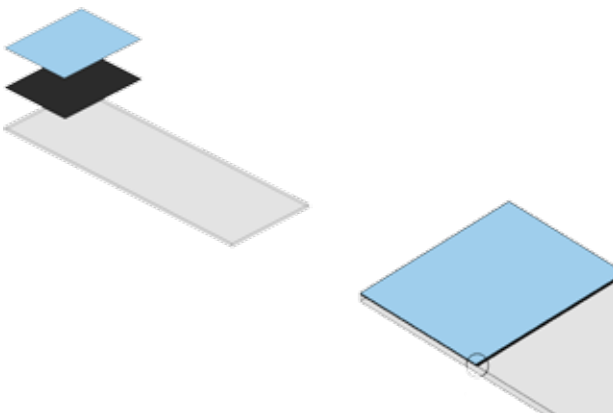






SHURMARK® MICROSCOPE SLIDES FOR PRIMERA SLIDE PRINTER



Different than the Thermo SlideMate™, the Signature® Slide Printer utilizes a specialty ribbon formulation, print head and slide surface to transfer its content. The imprints are resistant to reagent exposure and physical wear. The Signature® Slide Printer can print on the SHURMark® Routine Colored Microscope Slides and SHURMark® Contrasting Slides.

Warning: The Signature® Slide Printer will NOT accept 90 degree corner slides. Due to the printer's design, the 90 degree corner slides will jam when exiting the dispensing hopper.

SHURMark® Routine Colored Microscope Slides for use with Pens & Pencils					
Corners	 45° Clipped Corners		 Beveled Edges		
Colors	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	
 WHITE	SM-CS-45W	SM-CCS-45W	SM-CS-BW	SM-CCS-BW	
 BLUE	SM-CS-45B	SM-CCS-45B	SM-CS-BB	SM-CCS-BB	
 YELLOW	SM-CS-45Y	SM-CCS-45Y	SM-CS-BY	SM-CCS-BY	
 GREEN	SM-CS-45G	SM-CCS-45G	SM-CS-BG	SM-CCS-BG	
 LILAC	SM-CS-45L	SM-CCS-45L	SM-CS-BL	SM-CCS-BL	
 PINK	SM-CS-45P	SM-CCS-45P	SM-CS-BP	SM-CCS-BP	

SHURMark® Contrasting Black Color Under Primary Surface Color for use with Etching Labels					
Corners	 45° Clipped Corners		 Beveled Edges		
Colors	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	
 WHITE	SM-CU-45W	SM-CCU-45W	SM-CU-BW	SM-CCU-BW	
 BLUE	SM-CU-45B	SM-CCU-45B	SM-CU-BB	SM-CCU-BB	
 YELLOW	SM-CU-45Y	SM-CCU-45Y	SM-CU-BY	SM-CCU-BY	
 GREEN	SM-CU-45G	SM-CCU-45G	SM-CU-BG	SM-CCU-BG	
 LILAC	SM-CU-45L	SM-CCU-45L	SM-CU-BL	SM-CCU-BL	
 PINK	SM-CU-45P	SM-CCU-45P	SM-CU-BP	SM-CCU-BP	

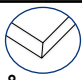
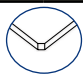







SHURMARK® MICROSCOPE SLIDES FOR LEICA/SAKURA PRINTERS

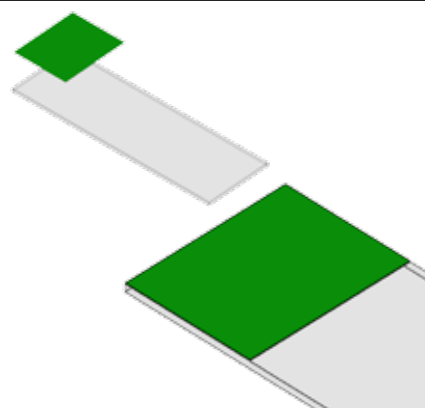


Both the Leica IP C and Sakura Tissue-Tek® AutoWrite® are ink printers. They utilize a specialty ink that produces imprints resistant to reagent exposure and physical wear. Ink printers can only print on the SHURMark® PLUS Color-Frosted Slides.

Warning: Never substitute with other slide types. The ink will not absorb properly, which will result in losing printed content after reagent exposure.


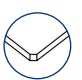


SHURMark® PLUS Color-Frosted Slides for use with Ink-Jet Labelers, Pens & Pencils

Corners		 90° Corners		 45° Clipped Corners		 Beveled Edges	
Colors		UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #
	WHITE	SMP-STT-W	SMP-SCTT-W	SMP-STT-45W	SMP-SCTT-45W	SMP-STT-BW	SMP-SCTT-BW
	BLUE	SMP-STT-B	SMP-SCTT-B	SMP-STT-45B	SMP-SCTT-45B	SMP-STT-BB	SMP-SCTT-BB
	YELLOW	SMP-STT-Y	SMP-SCTT-Y	SMP-STT-45Y	SMP-SCTT-45Y	SMP-STT-BY	SMP-SCTT-BY
	GREEN	SMP-STT-G	SMP-SCTT-G	SMP-STT-45G	SMP-SCTT-45G	SMP-STT-BG	SMP-SCTT-BG
	LILAC	SMP-STT-L	SMP-SCTT-L	SMP-STT-45L	SMP-SCTT-45L	SMP-STT-BL	SMP-SCTT-BL
	PINK	SMP-STT-P	SMP-SCTT-P	SMP-STT-45P	SMP-SCTT-45P	SMP-STT-BP	SMP-SCTT-BP

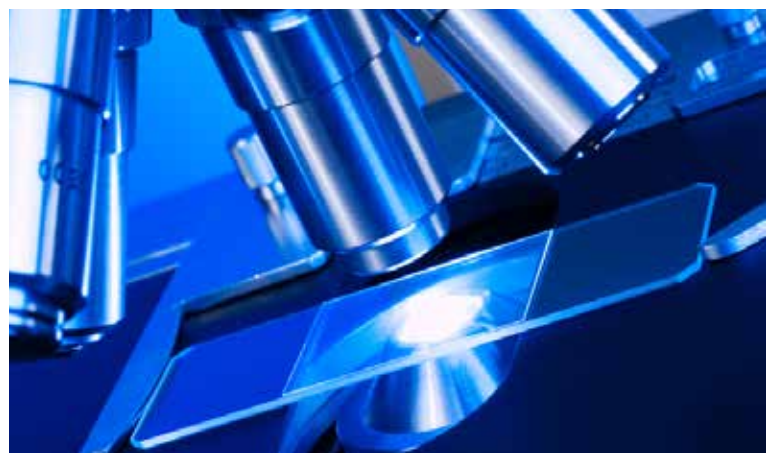


SHURMARK® MICROSCOPE PLAIN SLIDES

SHURMark® Frosted & Plain Microscope Slides

Corners	 90° Corners	 45° Clipped Corners	 Beveled Edges	
Style	UN-CHARGED Cat #	UNCHARGED Cat #	UNCHARGED Cat #	
Plain	SM-PTT	SM-PTT-45	SM-PTT-B	

SHURMark®
Beveled Edge
Slides are more
comfortable and
easier
to handle




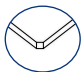

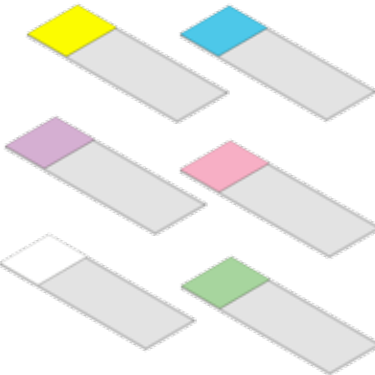






SHURMARK® MICROSCOPE SLIDES FOR THERMO SLIDEMATE PRINTERS






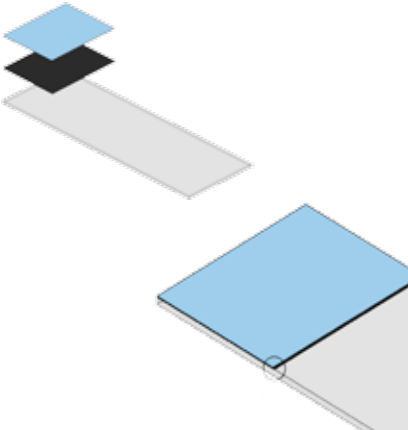






The Thermo SlideMate™ utilizes a specialty ribbon formulation, print head and slide surface to transfer its content. The imprints are resistant to reagent exposure and physical wear. The SlideMate™ can print on the SHURMark® Routine Colored, SHURMark® Contrasting and SHURMark® Frosted Slides.

Warning: SHURMark® PLUS Color-Frosted Slides cannot be substituted. When printed, it produces inadequate printing not resistant to chemical exposure or physical wear.




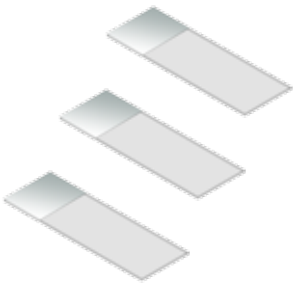
SHURMark® Routine Colored Microscope Slides for use with Pens & Pencils

Corners		 90° Corners		 45° Clipped Corners		 Beveled Edges		
Colors		UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	
	WHITE	SM-CS-W	SM-CCS-W	SM-CS-45W	SM-CCS-45W	SM-CS-BW	SM-CCS-BW	
	BLUE	SM-CS-B	SM-CCS-B	SM-CS-45B	SM-CCS-45B	SM-CS-BB	SM-CCS-BB	
	YELLOW	SM-CS-Y	SM-CCS-Y	SM-CS-45Y	SM-CCS-45Y	SM-CS-BY	SM-CCS-BY	
	GREEN	SM-CS-G	SM-CCS-G	SM-CS-45G	SM-CCS-45G	SM-CS-BG	SM-CCS-BG	
	LILAC	SM-CS-L	SM-CCS-L	SM-CS-45L	SM-CCS-45L	SM-CS-BL	SM-CCS-BL	
	PINK	SM-CS-P	SM-CCS-P	SM-CS-45P	SM-CCS-45P	SM-CS-BP	SM-CCS-BP	

SHURMark® Contrasting Black Color Under Primary Surface Color for use with Etching Labelers

Corners		 90° Corners		 45° Clipped Corners		 Beveled Edges		
Colors		UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	
	WHITE	SM-CU-W	SM-CCU-W	SM-CU-45W	SM-CCU-45W	SM-CU-BW	SM-CCU-BW	
	BLUE	SM-CU-B	SM-CCU-B	SM-CU-45B	SM-CCU-45B	SM-CU-BB	SM-CCU-BB	
	YELLOW	SM-CU-Y	SM-CCU-Y	SM-CU-45Y	SM-CCU-45Y	SM-CU-BY	SM-CCU-BY	
	GREEN	SM-CU-G	SM-CCU-G	SM-CU-45G	SM-CCU-45G	SM-CU-BG	SM-CCU-BG	
	LILAC	SM-CU-L	SM-CCU-L	SM-CU-45L	SM-CCU-45L	SM-CU-BL	SM-CCU-BL	
	PINK	SM-CU-P	SM-CCU-P	SM-CU-45P	SM-CCU-45P	SM-CU-BP	SM-CCU-BP	

SHURMark® Frosted Microscope Slides

Corners		 90° Corners		 45° Clipped Corners		 Beveled Edges		
Style		UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	UNCHARGED Cat #	CHARGED Cat #	
Frosted One End, One Side		SM-SSF	SM-SSFC	SM-SSF-45	SM-SSFC-45	SM-SSF-B	SM-SSFC-B	
Frosted One End, Two Sides		SM-DSF	SM-DSFC	SM-DSF-45	SM-DSFC-45	SM-DSF-B	SM-DSFC-B	

PLASTIC SLIDE BOXES WITH HINGED COVERS



Slide Boxes - Hinged Cover

10 Position Slide Box, ABS.
2pcs/ea; 10ea/pk, 5pk/cs

Cat. #	Color
SB-10P-O	Orange
SB-10P-G	Green
SB-10P-B	Blue
SB-10P-W	White



Slide Boxes - Hinged Cover

25 Position Slide Box w/ Cork Insert, ABS.
2pcs/ea; 10ea/pk, 5pk/cs

Cat. #	Color
SB-25PC-O	Orange
SB-25PC-G	Green
SB-25PC-B	Blue
SB-25PC-W	White



Slide Boxes - Hinged Cover

50 Position Slide Box, ABS.
1pcs/ea; 5ea/pk, 25pk/cs

Cat. #	Color
SB-50P-OH	Orange
SB-50P-GH	Green
SB-50P-BH	Blue
SB-50P-WH	White



Slide Boxes - Hinged Cover

100 Position Slide Box, ABS.
1pcs/ea; 5ea/pk, 10pk/cs

Cat. #	Color
SB-100P-O	Orange
SB-100P-G	Green
SB-100P-B	Blue
SB-100P-W	White



Slide Boxes - w/Metal Clasp

100 Position Slide Box, ABS.

Cat. #	Color
SB-100-O-M	Orange
SB-100-G-M	Green
SB-100-B-M	Blue
SB-100-W-M	White

Wooden Slide Boxes with Hinged Covers



Slide Boxes - with Metal Hinge & Clasp

1pcs/ea; 5ea/pk, 10pk/cs

Cat. #	Slide Positions
SB-50-W	50
SB-100-W	100



Plastic Slide Boxes with Removable Lids

Slide Boxes - Removable Lids

25-Position Slide Boxes, ABS.
50-Position Slide Boxes, ABS.
2pcs/ea; 10ea/pk, 5pk/cs

Cat. #s		
	50- Position	Color
SB-25P-O	SB-50P-O-S	Orange
SB-25P-G	SB-50P-G-S	Green
SB-25P-B	SB-50P-B-S	Blue
SB-25P-W	SB-50P-W-S	White



PLASTIC SLIDE MAILERS

SHUR**Mark**® Plastic Slide Mailers are made of shock-resistant ABS for secure shipping and storage. Each mailer accepts 25 x 75mm or 1" x 3" microscope slides; 10ea/pk, 20pk/cs



SM-1P



SM-2P



SM-5P

Cardboard Slide Mailers & Folders

SHUR**Mark**® High-quality Cardboard Slide Mailers & Folders are made of heavy-duty white cardboard. Each position features a finger-slot for easy slide removal and is sized to accept 25 x 75mm or 1 x 3" microscope slides. 10, 20, & 30 position slide folders have indexes for recording slide data; 10ea/pk, 5pk/cs



SM-1C



SM-2C



SM-3C



SM-4C



SM-5C



SH-10C



SH-20C

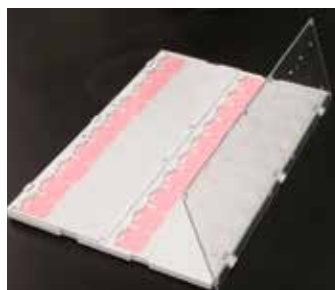


SH-30C

Plastic Slide Holders & Trays

ClearView Slide Holder

20-Position Slide Holder with double-hinged clear cover, ABS. 5ea/pk, 4pk/cs



Cat. #	Color
CVSH-W	White



Slide Trays

20-Position Tray, ABS. 5ea/pk, 4pk/cs

Cat. #	Color
SF-20P-Y	Yellow
SF-20P-G	Green
SF-20P-B	Blue
SF-20P-W	White

what's NEW?



Customer Focus

- Customer Solutions
- Technical Support

NEW Company

Innovation

- NEW Products
- NEW Solutions
- Proven Technology

NEW Solutions

Continuous Improvement

- Specimen & Employee Safety
- Quality Control
- Certifications

Reliability

- Throughput
- Dependability
- Accountability

NEW Products

Microscope Slide Staining Dishes & Accessories

Glass Staining Dishes

Cat. #	Description
SD-5-G	Vertical Staining Dish, Glass, 5-positions; 1pcs/ea; 5ea/pk, 10pk/cs
SD-10-G	Horizontal Staining Dish, Glass, 10-positions; 1pcs/ea; 5ea/pk, 10pk/cs



Plastic Staining Dishes & Rack

Cat. #	Description
SS-WLW	Staining Dish, 24-position, ABS, WHITE; 15/cs
SS-WLG	Staining Dish, 24-position, ABS, GREEN; solvent resistant; 15/cs
SS-SRG-24	Staining Slide Rack, 24-position, ABS, GREY; 5/pk



Microscope Slide Apparatus



Cat. #	Description
SDS	Slide Drying Stand, 77 slide capacity, white, HDPE, stain resistant. (0.5 x 7 x 11.75) in.; 1pcs/ea
MSD	Microscope Slide Dispenser (holds up to 72 Slides); 1/ea
SM-SG	SHUR/Mark® Slide Grip (holds up to 5 slides), chemical and solvent resistant, capable of withstanding temps up to 80°C; 6ea/cs

SHURMark® Marking Pens

For cassettes and slides, solvent resistant, fast drying, won't smear. The only laboratory marking pen with priming action and replaceable felt tip.

Cat. #	Description
MP-B	SHURMark® Marking pen, BLACK, for labeling slides or cassettes; 12ea/pk, 5pk/cs
MP-R	SHURMark® Marking pen, RED, for labeling slides or cassettes; 12ea/pk, 5pk/cs
MP-RT	SHURMark® Replacement Tips (for marking pens); 12ea/pk, 5pk/cs



SHURStor™ Slide Storage File

Convenient, space saving and portable. Stores 20,000 slides or 5,000 cassettes/case. Sturdy, rugged, moisture resistant and is stackable up to 10 files high.

Cat. #	Description
SF-S	SHURStor™ Slide Storage File w/4 drawers, capacity 2000 slides/file; 10 stackable files/cs
SF-B	SHURStor™ Tissue Block Storage w/8 drawers, capacity 500 cassettes/file; 10 stackable files/cs



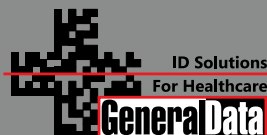
More Information

PH: 844.643.1129

www.general-data.com/hc

General Data Healthcare helps labs enhance productivity and improve workflow with innovative products, solutions and support that focus on delivering better patient safety & care.

©2014 General Data Healthcare, Inc. All rights reserved. Specifications are current at the time of publication; however, they are subject to change without notification. All trademarks and registered trademarks are property of General Data Healthcare, Inc., Cincinnati, Ohio, USA.



General Data Healthcare, Inc.
Headquarters: 4354 Ferguson Drive, Cincinnati, OH 45245
TBS Facility: 3014 Croasdaile Drive, Durham, NC 27705

Sales Inquiries:
PH: 844.643.1129
Email: talktous@general-data.com
Web: www.general-data.com/hc

Service:
PH: 800.245.5029
Email: customersupport@general-data.com