

*Histology Innovation for a NEW Generation*

**Staining**

# SHURStain™ Manual Stainers



**Proven Technology. Superior Reliability.  
Protecting Every Life Story In Your Lab.**

Each tissue specimen your lab processes has a life story behind it. Your mission is to ensure all your specimens receive optimal processing in order to deliver the best possible care. SHURStain™ Manual Stainers enable your lab to deliver increased productivity, advanced specimen safety and highly accurate results.

## Customizable Manual Stainers Deliver Enhanced Throughput and Accuracy

There's a patient's life story behind every tissue specimen your lab processes, and it is your job to ensure every life story receives the best possible care. The SHURStain™ utilizes customizable manual stainer racks that are designed for the needs of today's lab - where increased reliability, efficiency, throughput, accuracy and safety are all critical factors for delivering the highest quality care to your patients.

### ✓ Preferred Choice for Smaller Volume Labs


These modular stainers enable you to select the number of individual staining positions necessary to meet your histology, cytology, dermatology, immunohistochemical, hematology, microbiology or special staining application. They are produced in 3, 6, & 12 well configurations, and multiple units can be interlocked back-to-back to increase staining volume capabilities.

### ✓ Easy Access With Minimal Solvent Loss

Sliding lids keep stains and reagents covered when not in active use, preventing unwanted evaporation.


### ✓ Interlocking Design Allow Multiple Units to be Attached Together

Shown below is the SS-6L interlocked back-to-back.

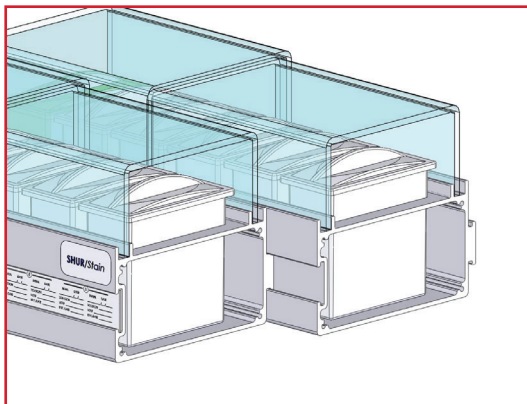

## SS-12L

12 position  
Manual Stain Rack

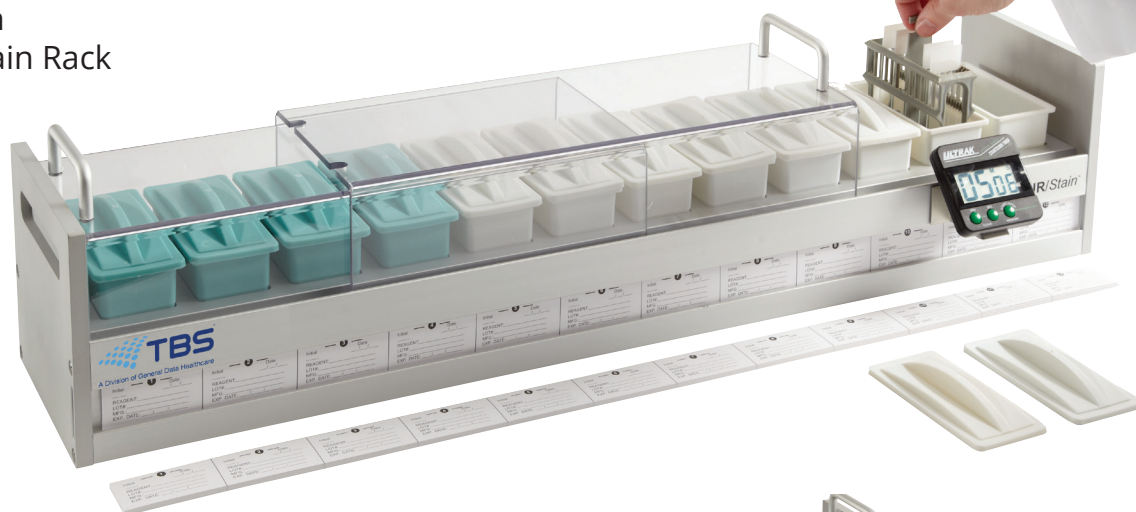

## SS-3L

3 position  
Manual Stain Rack




## SS-6L

6 position  
Manual Stain Rack



## Stainer Accessories

1. Digital Timer attaches directly to the stainer with convenient bracket.
2. Green & White 250ml Replacement Wells; shown with 24-position Solvent Resistant Slide Carrier
3. Built-In cut-out handles to conveniently move stainers around the lab.
4. Sliding Lids opens from either end to minimize reagent evaporation. No need to replace individual lids during daily use.

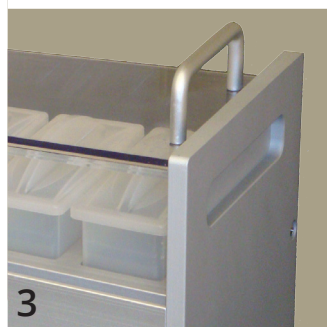

1		2		3	
Initial	Date	Initial	Date	Initial	Date
AC	07/14/2009				
REAGENT:	HEM	REAGENT:		REAGENT:	
LOT#:	29007A	LOT#:		LOT#:	
MFG.:	TBS	MFG.:		MFG.:	
EXP. DATE:	Feb 6, 2011	EXP. DATE:		EXP. DATE:	

Each SHURStain™ rack comes with a supply of replaceable documentation strips to record reagent brand, lot #, date changed, etc. The number of data positions is proportional to the number of wells.

# Product Specifications:

Manual Stainer Rack	SS-3L	SS-6L	SS-12L
Dimensions: DxWxH (in)	(6x8x7)(160x203-x178)(in)(mm)	(6x15x6)(160x375x156)(in)(mm)	(6x28x6)(160x718x156)(in)(mm)

## How To Order:

SS-3L	SHURStain™ Manual Stain Rack, 3 position; includes 2 white 250ml wells (SS-WLW), 1 green solvent resistant 250 ml well (SS-WLG), 1ea. 24-position plastic slide carrier (SS-SRG-24). Individual lids are provided for use during extended inactivity & a see-through lexan cover is included for convenience & safety during periods of frequent use. Multiple racks can be attached back to back. Digital timer and (12) label strips included (SS-LS-3L ) to record reagent expiration dates, lot #'s, etc.
SS-6L	SHURStain™ Manual Stain Rack, 6 position; includes 4 white 250ml wells (SS-WLW), 2 green solvent resistant 250 ml well (SS-WLG), 1ea. 24-position plastic slide carrier (SS-SRG-24). Individual lids are provided for use during extended inactivity & see-through lexan cover is included for convenience & safety during periods of frequent use. Multiple racks can be attached back to back. Digital timer and (12) label strips included (SS-LS-6L) to record reagent expiration dates, lot #'s, etc.
SS-12L	SHURStain™ Manual Stain Rack, 12 position; includes 8 white 250ml wells (SS-WLW), 4 green solvent resistant 250 ml well (SS-WLG), 1ea. 24-position plastic slide carrier (SS-SRG-24). Individual lids are provided for use during extended inactivity & see-through lexan cover is included for convenience & safety during periods of frequent use. Multiple racks can be attached back to back. Digital timer and (12) label strips included (SS-LS-12L) to record reagent expiration dates, lot #'s, etc.
SS-SRG-24	Slide Rack, plastic, Grey, 24-position, solvent resistant
SS-WLW	Well, 250ml, w/ lid, white
SS-WLG	Well, 250ml. w/ lid, green, xylene resistant
SS-DT	Digital Timer, countdown, built-in memory, loud alarm, large display
SS-LS-3L	Label Strip, replacement strip for SS-3L with 3 data entry positions, 8", 12 ea/pk
SS-LS-6L	Label Strip, replacement strip for SS-6L with 6 data entry positions, 15", 12 ea/pk
SS-LS-12L	Label Strip, replacement strip for SS-12L with 12 data entry positions, 28", 12 ea/pk

## Related Products:

Instrumentation:	SHURDry™ Slider Dryers II & III, SHURMark® Slide Labeler
Consumables:	SHURMark® Microscope Slides & Cover Glass, Microscope Slide Drying Stand, SHURMark® Microscope Slide Folders, Boxes & Mailers, SHURStor™ Slide Storage Files
Reagents:	SHURStain™ Routine & Special Stain Kits


### More Information

PH: 844.643.1129

[www.general-data.com/hc](http://www.general-data.com/hc)

*General Data Healthcare helps labs enhance productivity and improve workflow with innovative products, solutions and support that focus on delivering better patient safety & care.*

©2014 General Data Healthcare, Inc. All rights reserved. Specifications are current at the time of publication; however, they are subject to change without notification. All trademarks and registered trademarks are property of General Data Healthcare, Inc., Cincinnati, Ohio, USA.

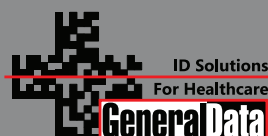

General Data Healthcare, Inc.  
Headquarters: 4354 Ferguson Drive, Cincinnati, OH 45245  
TBS Facility: 3014 Croasdaile Drive, Durham, NC 27705

Sales Inquiries:  
PH: 844.643.1129  
Email: [talktous@general-data.com](mailto:talktous@general-data.com)  
Web: [www.general-data.com/hc](http://www.general-data.com/hc)

Service:  
PH: 800.245.5029  
Email: [customersupport@general-data.com](mailto:customersupport@general-data.com)